





Sterling Couch Range

Therapy

Couch Range



Patient Trolley Range



ISO 13485

- 2016 -

Seating Range





SEERS MEDICAL Ltd Kenton Road, Debenham, Suffolk IP14 6LA, United Kingdom © Copyright SEERS Medical Ltd 2022



Tel +44 (0)1728 861 590 **Fax** +44 (0)1728 861 868

ISO 14001

- 2015 -

ISO 9001

2008

info@seersmedical.com www.seersmedical.com

> lssue № 8/22 Part № 06Z-8005-SF

Company Reg. № 6572698 VAT № 931596702 Printed in the United Kingdom







Clinical Range

Clinical Couches for Patient Examination and Minor Treatment Procedures




Introduction

In this new brochure we cover our wide range of patient examination furniture for various medical procedures such as Gynaecology, Ultrasound, Phlebotomy, Lymphoedema and Podiatry. There is a wide choice of models and accessories to cater for various examination procedures as well as specialist treatment and minor operation techniques. Most models are variable in height with a choice of both manual or electric positioning and height operation. The development of these products has been achieved by applying our many years of experience and knowledge in this market. Much thought has been given to the design, to benefit both the comfort and security of the patients as well as providing manual handling solutions which promote good posture and ergonomics to all clinical staff when they are using this equipment. With the management and control of our own assembly line at our premises in Debenham, Suffolk in the UK, our couches comply to the CE medical devices directive 93/42/EEC, ensuring we supply high quality products every time, with short lead times to meet the tightest of deadlines. We trust you will experience this yourself when you purchase a **SEERS** Medical product. **SEERS Medical are ISO9001 and ISO14001 accredited**.

Contents

4-5	Medicare 2 Section Couch
6-7	Medicare 3 Section Couch
8-9	Medicare Mobile 2 Section Couch
10-11	Medicare Mobile 3 Section Couch
12-13	Innovation Deluxe 2 Section Couch
14-15	Innovation Deluxe 3 Section Couch
16-17	Innovation Deluxe Daycare Couch
18-19	Innovation Deluxe Dialysis Couch
20-21	Medicare Multi-Couch - Single Footr
22-23	Medicare Podiatry Couch
24-25	Medicare Orthopaedic Couch
26-27	Medicare Phlebotomy Couch
28-29	Medicare Echocardiography Couch
30-31	Medicare Ultrasound Couch
32-33	Medicare Gynaecology Couch
34-35	Medicare Deluxe Gynaecology Couc
36-37	Medicare GP Gynaecology Couch
38-39	Medicare Bariatric 2 Section Couch
40-41	Medicare Bariatric 3 Section Couch
42-43	Medicare Bariatric Drop End Couch
44-45	Shower Trolley
46-47	Optional Accessories
48-49	Chairs, Stools & Vinyl Colour Options

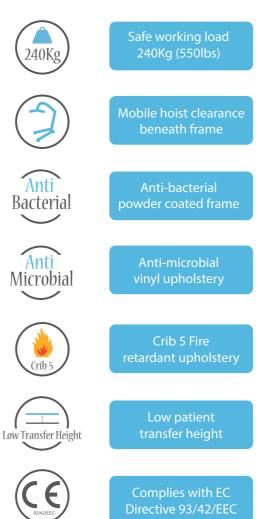

Technical Information

At **SEERS** Medical we are using state of the art design tools, such as 3D solid modelling and in house rapid prototyping, to assist us in getting new designs into development and manufacture. Whilst designing new products, all manufacturing capabilities available are taken into consideration, such as CNC wood cutting, laser profiling and robot welding, which deliver reliable build quality and repeatability of accurate components, time after time. A typical innovative use of these processes is evident in a high loaded components such as the ram bracket, a common item on many **SEERS** couches, where using laser profiling we have been able to make the parts interlock, adding immense strength, even before welding! This is typical of the attention to detail that we pride ourselves with, and sets us apart from other couch manufacturers, giving us confidence in the reliability and longevity of the products we manufacture.

Mechanical Safety & Load Testing

All **SEERS** Medical products have been tested to the applicable EN 60601-1:2006+A12:2014 (E) the General requirements for basic safety and essential performance. Our couches have been thoroughly tested with a two and a half times safe working load (in excess of 600Kg loaded evenly, for in excess of one hour, to prove load support capability without failure. Our products are tested in accordance to Clauses 9.8.3.2-3 Static Forces of EN60601-1:2006 Not only do our products have high safe working loads (typically 240Kg), but our products also lift the rated safe working load from the minimum height to the maximum height smoothly.

Electrical Safety

All **SEERS** Medical products comply to EN 60601-1:2006+A12:2014 (E) the General requirements for basic safety and essential performance. In addition, the actuation and control systems used on our products are marked with the UL mark, demonstrating compliance to UL60601-1 and CAN/CSA-C22.2 No. 601.1-M90, meeting the requirements for USA and Canada. The actuation control system is also marked with approval to C-ETL-US No. 3154358.

CE Marking & Quality System

SEERS Medical operates a Quality Management System in accredited and approved to ISO 13485:2016. **SEERS** Medical products comply with the Essential Requirements in Annex I of EC Directive 93/42/EC, with all relevant design, production records, risk assessments, product testing and labelling information recorded in the documentation file labelled the "Technical File". All products carry a unique serial number for trace-ability of build date, model number and associated build records.

Flammability

The upholstered couch sections satisfies flammability requirements to the standards of BS 7176:1995 (medium hazard) using ignition sources 0 (cigarette), 1 (gas) as defined in BS EN. 1921-1 and 2 and 5 (crib) BS 5852:1990.


Warranty

The variable height couches detailed in this brochure have a lifetime frame warranty and five years warranty on all other components.





CAD design technology

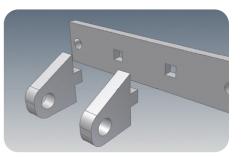


Interlocking structural components



Laser profiling



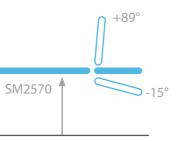



Height range 45cm to 98cm

Medicare 2 Section Couch

Anti Bacterial 240Kg Low Transfer He

A very strong and rigid couch available with a wide range of accessories. Suitable for supporting very heavy patients for examination and minor treatment procedures.



Standard Features

- Safe working load and lifting capacity of 240Kg (528lb
- Gas assisted backrest adjustable from -25° to 75° on
- Choice of hydraulic or electric height elevation
- Retractable wheel design with adjustable foot for stab
- Clearance beneath frame for mobile hoist
- Hand switch fitted as standard on electric models
- Frames pre-drilled for addition of accessories
- 63cm wide Easy Clean Upholstery as standard
- Choice of Standard or Standard+ upholstery colours




Low patient transfer height

STANDARD STANDARD+ **Available** Upholstery **Colours:** Black White ——— Light Grey ——— Cream —— Sky Blue Canard Dark Grev Olive Dark Blue — Lotus Green — Pepper

	Opt	ional Accessories
s)	9002	Side support rails (Pair)
	9002P	Side support cushions (Pair)
n	6004	Breathing hole and plug
0	6005	Paper roll holder
	6006	Upper limb support
	6007	Removable head/knee support
ility	6053	Protective base frame cover
	6054	Telescopic IV pole
	6057	Battery back-up facility
	-EPFS	Perimeter Foot Switch for SM2560
	-1WF	Single foot switch for SM2560
	-2WF	Double foot switch for SM2570
	9200/70	70cm Width Easy Clean Upholstery



Negative head section to -25°

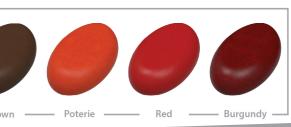







Medicare 3 Section Couch

Anti Bacterial 240Kg Low Transfer Hei

The Medicare 3 Section Couch provides a very stable and versatile support surface for prone, supine and seated treatment procedures and examinations. Both the elevating sections can be set to a wide range of angles for specific patient positioning.

Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)

- Gas assisted sections incorporate anti-trap design
- Choice of hydraulic or electric height elevation

- Purpose made precision bearings on all pivot points



Available Upholstery **Colours:**



6

SM3550

SM3560

SM3570


Standard Dimensions


*When fitted with 9300/70 Easy Clean Upholstery

Optional Acces	CORIDE

9003	Side support rails (Pair)
9002P	Side support cushions (Pair)
6004	Breathing hole and plug
6005	Paper roll holder
6006	Upper limb support
6007	Removable head/knee support
6053	Protective base frame cover
6054	Telescopic IV pole
6057	Battery back-up facility
-EPFS	Perimeter Foot Switch for SM3560
-1WF	Single foot switch for SM3560
-2WF	Double foot switch for SM3570
9300/70	70cm Width Easy Clean Upholstery



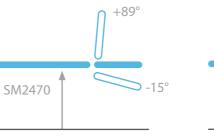

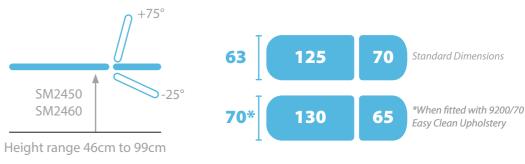

Height range 46cm to 99cm

Medicare Mobile 2 Section Couch

Anti Bacterial 240Kg Low Transfer He

Very strong and rigid 2 section mobile treatment couches supplied with side support rails and a paper roll holder. The large wheel system incorporates a steering facility and the couches have excellent height ranges with a low patient transfer height of just 46cm.

Models

Independently Braking Castors (IBC)

SM2450	Hydraulic with gas assisted backrest
SM2460	Electric with gas assisted backrest
SM2470	Electric with electric backrest
Central Lo	cking Castors (LMWD)
SM2650	Hydraulic with gas assisted backrest
SM2660	Electric with gas assisted backrest
SM2670	Electric with electric backrest



Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)
- Gas assisted backrest adjustable from -25° to 75° on SM2450, SM2460, SM2650 and SM2660
- Gas assisted backrest incorporates anti-trap design
- Choice of hydraulic or electric height elevation
- Excellent height range from 46cm to 99cm
- Side Rails & Paper Roll Holder fitted as standard
- Choice of independently or centrally locking castors
- Clearance beneath frame for mobile hoist
- Hand switch fitted as standard on electric models
- Frames pre-drilled for addition of accessories
- Anti-bacterial powder coated frame
- 63cm wide Easy Clean Upholstery as standard
- Choice of Standard or Standard+ upholstery colours



Model SM2650 Shown with optional 6053 Base Cover, 9002P Side Cushions, Colour Sky Blue

Available Upholstery **Colours:**




Optional Accessories

9002P Side support cushions (Pair) **6004** Breathing hole and plug Upper limb support 6006 6007 Removable head/knee support **6053** Protective base frame cover **6054** Telescopic IV pole **6057** Battery back-up facility -**1WF** Single foot switch for SM2X60 -**2WF** Double foot switch for SM2X70 **9200/70** 70cm Width Easy Clean Upholstery



Seated position 6006 Upper Limb Support



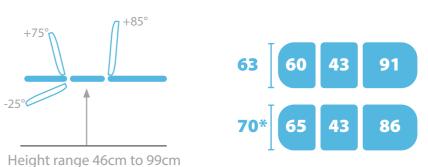



Medicare Mobile 3 Section Couch

Electric with gas assisted backrest

Hydraulic with gas assisted backrest

Electric with gas assisted backrest

Electric with electric backrest



Very strong and rigid 3 section mobile treatment couches supplied with side support rails and a paper roll holder as standard. The Large Mobile Wheel Design incorporates central locking castors & steering facility, and the couches have excellent height ranges with a low patient transfer height.

Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)

- Choice of hydraulic or electric height elevation
- Choice of independently or centrally locking castors
- Side Rails & Paper Roll Holder fitted as standard
- Purpose made precision bearings on all pivot points
- Clearance beneath frame for mobile hoist

- Anti-bacterial powder coated frame
- 63cm wide Easy Clean Upholstery as standard
- Choice of Standard or Standard+ upholstery colours

Model SM3650 Colour Sky Blue



240Kg

Models

SM3450 SM3460

SM3470

SM3650

SM3660

SM3670

Independently Braking Castors (IBC)

Central Locking Castors (LMWD)

Available

Colours:

Standard Dimensions


*When fitted with 9300/70 Easy Clean Upholstery

Optional Accessories


9002P	Side support cushions (Pair)
6004	Breathing hole and plug
6006	Upper limb support
6007	Removable head/knee support
6053	Protective base frame cover
6054	Telescopic IV pole
6057	Battery back-up facility
-1WF	Single foot switch for SM3X60
-2WF	Double foot switch for SM3X70
9300/70	70cm Width Easy Clean Upholstery



Low patient transfer height

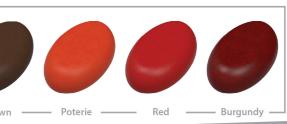



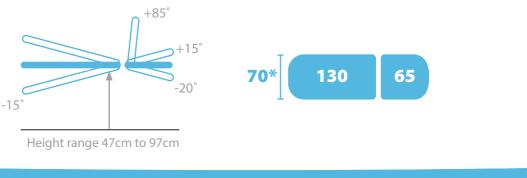

Innovation Deluxe 2 Section



The Innovation Deluxe 2 section couch is suitable for use within specialist critical care environments, endoscopy, ultrasound scanning and minor treatment procedures. Electric height, backrest and 2-way tilt ensure smooth transition when positioning patients for specialist procedures and treatments. Auto-electric CPR functionality and auto-levelling software return the couch quickly to a fully horizontal resuscitation position in the event of an emergency.

Models

MG2475	Electric with electric backrest and 2-w
	tilt. Large independently locking
	wheel design
MG2675	Electric with electric backrest and 2-w
	tilt. Large central locking wheel desigr
	with steering facility



Standard Features

- Electric trendelenburg tilt to -/+15°
- - Safe working load and lifting capacity of 250Kg (550lbs)
 - Auto-electric CPR function via handset
 - Auto-intelligent operation
 - Vertical lifting design to minimise space and provide maximum stability throughout all profiling positions
 - Smooth moulded base cover minimises dirt traps and

 - Large TENTE® Ø125mm castors
 - **SEERS** ComfortPlus[™] Upholstery features anti-microbial vinyl incorporating pressure relieving design for additional patient comfort with a foam depth of 60mm
- 70cm wide upholstery
- Clearance beneath frame for mobile hoist
- Head Cushion & Paper Roll Holder fitted as standard
- Choice of Standard or Standard+ upholstery colours

Model MG2675 with large central locking wheel design with steering facility





Negative head section to -20°

Available Upholstery **Colours:**



Optional Accessories

6003	Side support rails (Pair)
9002P	Side support rail cushions (Pair
6004	Breathing hole and plug
6006	Upper limb support
6054	Telescopic IV pole
6057	Battery back-up facility



Trendelenburg tilt to -15°





Innovation Deluxe 3 Section



The Innovation Deluxe 3 section couch is designed for use within specialist critical care environments. Electric height, backrest and foot section profiling and 2-way tilt ensure a smooth transition when positioning patients for specialist procedures and treatments. Auto-electric CPR functionality and auto levelling software return the couch quickly to a fully horizontal resuscitation position in the event of an emergency.

Models

Μ	-	л	7	
111	 61			
	-			-

MG3675

foot section with auto-levelling, 2-way tilt and CPR facility. Large independently Electric profiling of height, backrest and tilt and CPR facility. Large central locking wheel design with steering facility

Electric profiling of height, backrest and

Standard Features

- Electric trendelenburg tilt to -/+15°

- Safe working load and lifting capacity of 250Kg (550lbs)
- Auto-electric CPR function via handset
- Auto-intelligent operation
- Vertical lifting design to minimise space and provide maximum stability throughout all profiling positions
- Smooth moulded base cover minimises dirt traps and facilitates easy cleaning
- - Large TENTE[®] Ø125mm castors
 - SEERS ComfortPlus[™] Upholstery features anti-microbial vinyl incorporating pressure relieving design for additional patient comfort with a foam depth of 60mm

 - Clearance beneath frame for mobile hoist
 - Head Cushion & Paper Roll Holder fitted as standard
 - Choice of Standard or Standard+ upholstery colours


Model MG3475 with large independently braked wheel design




Low patient transfer height

Available Upholstery **Colours:**



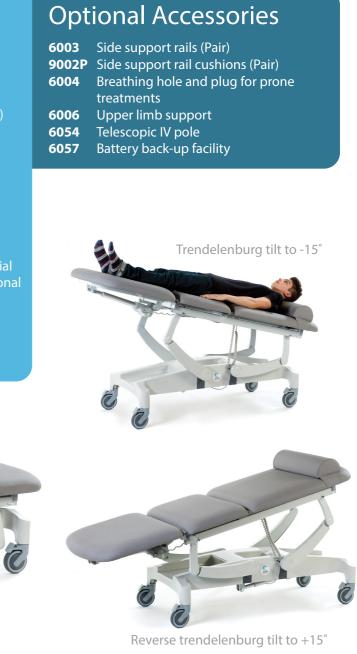



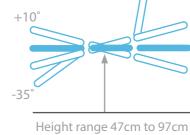

Innovation Deluxe Daycare



The Innovation Deluxe Daycare couch is designed for specialist medical procedures, critical care and oncology units. The fully electric functionality provides smooth lifting and patient profiling during patient examination and treatment. The couch features 2-way tilt and auto-electric CPR facility with auto-levelling software.

Models

MG3493

MG3693

Electric profiling of height, backrest and foot section with auto-levelling, 2-way tilt and CPR facility. Large independently Electric profiling of height, backrest and tilt and CPR facility. Large central locking wheel design with steering facility

Sky Blue -

Canard _____

— Dark Blue -



Dark Grev

Black

Standard Features

- Electric trendelenburg tilt to -/+15°

- Safe working load and lifting capacity of 250Kg (550lbs)
- Auto-electric CPR function via handset
- Auto-intelligent operation
- Vertical lifting design to minimise space and provide maximum stability throughout all profiling positions
- Smooth moulded base cover minimises dirt traps and facilitates easy cleaning
- - Large TENTE[®] Ø125mm castors
 - **SEERS** ComfortPlus[™] Upholstery features anti-microbial vinyl incorporating pressure relieving design for additional patient comfort with a foam depth of 60mm


 - Clearance beneath frame for mobile hoist
 - Head Cushion, Paper Roll Holder, Side Rails and Side Support Cushions fitted as standard
 - Fold down side rails and cushions provide patient comfort and security



Colours:

2-way trendelenburg tilt

STANDARD STANDARD+ White — Light Grey — Cream — Olive —— Lotus Green —— Pepper

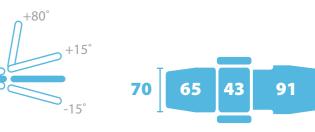

Optional Accessories

6057 6241

6054 Telescopic IV pole Battery back-up facility Patient lap tray 6242 Adjustable foot board



Model 6241 optional


patient lap tray




Model 6242 optional adjustable foot board



Low patient transfer height

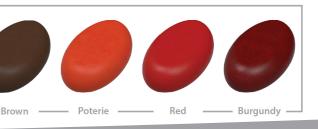



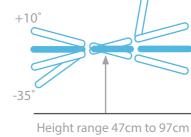

Innovation Deluxe Dialysis



The Innovation Deluxe Dialysis couch has been designed for use within specialist renal dialysis environments. The couch features auto-electric CPR and fully electric profiling with a 2-way tilt facility to ensure a smooth transition when positioning patients for dialysis procedures and treatments.

Models

MG3490

MG3690

Electric profiling of height, backrest and and CPR facility. Large independently braked wheel design Electric profiling of height, backrest and tilt and CPR facility. Large central locking wheel design with steering facility

> Model MG3490 with large independently braked wheel design

Standard Features

- Electric trendelenburg tilt to -/+15°

- Auto-electric CPR function via handset
- Vertical lifting design to minimise space and provide
- Smooth moulded base cover minimises dirt traps and facilitates easy cleaning
- - patient comfort with a foam depth of 60mm


 - board included for patient comfort




Fully profiled for patient comfort

Available Upholstery **Colours:**



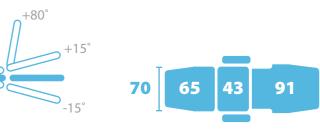



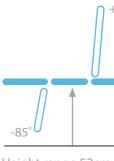

Medicare Multi-Couch

Anti Bacterial (**2**40Kg

Incorporating both an adjustable angle backrest and a single drop end foot section, these models are suitable for a wide range of seated and supine patient examination procedures. The slightly narrow foot section enables the addition of leg supports or stirrups for gynaecological procedures.

Models

SM9556 SM9566 SM9576

Electric with electric backrest



Model 6045 optional leg supports for gynaecology

Model SM9566 Colour Light Grey

Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)

- Gas assisted sections incorporate an anti-trap design
- Choice of hydraulic or electric height elevation

- Purpose made precision bearings on all pivot points
- Anti-bacterial powder coated frame
- - 65cm standard width upholstery
 - Choice of Standard and Standard+ upholstery colours





+80



\frown		
()ntio	nal Acc	essories

80°	

6004	Breathing hole and plug
6005	Paper roll holder
6006	Upper limb support
6007	Removable head/knee cushion
6044	Lithotomy stirrups
6045	Leg supports for gynaecology
6052	Protective base frame cover
6054	Telescopic IV pole
6057	Battery back-up facility
6057 -1WF -2WF	





Low patient transfer height



Model 6044 optional lithotomy stirrups



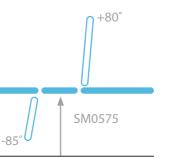

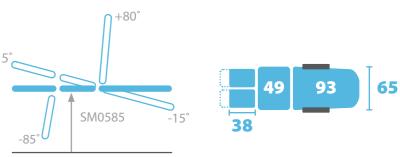

Standard Features

• Footrests extend for tall patients

• Safe working load and lifting capacity of 240Kg (528lbs)

incorporate anti-trap design gas struts

• Fold-down armrests have angular adjustment

Height range 52cm to 106cm

Medicare Podiatry Couch

Anti Bacterial (📕 240Kg

Incorporating both an adjustable angle backrest with patient armrests and dual adjustable footrests, these models are suitable for a wide range of seated examination procedures such as podiatry, leg ulcer clinics, orthopaedic casting (plaster room), as well as being ideal for GP surgeries and health centres.

Models

SM0575 SM0585

Electric height and backrest



6004 Breathing hole and plug for prone treatments



Available Upholstery Colours:



Optional Accessories

6052	Protective base frame cover
6007	Head Cushion (Included on SM0585)
6005	Paper roll holder (Included on SM0585)
6004	Breathing hole and plug
6054	Telescopic IV Pole Attachment
6057	Battery back-up facility
-2WF	Double foot switch for SM0575
-3WF	Triple foot switch for SM0585




Matching headrest as standard on Model SM0585. Option on Model SM0575



Low patient transfer height



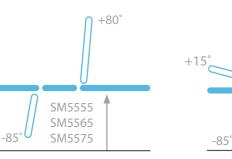

Height range 52cm to 106cm

Medicare Orthopaedic Couch

Designed specifically for plaster casting technicians, this couch features individually adjustable leg sections and is available with a range of optional accessories to assist and support the technician in carrying out various casting procedures. The

Low Transfer Height

Deluxe version is supplied as a package with base shroud, set of accessories and matching stool.

Models

SM5555	
SM5565	
SM5575	
SM5585	
SM5585D	

Hydraulic with gas assisted backrest Electric with electric backrest Electric with electric height, tilt and backrest Deluxe model with stool, base cover and Orthopaedic accessories set

Model SM5565 Colour Dark Blue Anti Bacterial

240Kg

Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)
- Gas assisted or electric backrest adjustable from 0° to +
- Trendelenburg tilt to -15° on Models SM5585/ SM5585D
- Gas assisted sections incorporate anti-trap design
- Footrests extend for tall patients
- Choice of hydraulic or electric height elevation
- Compact frame design with retractable wheel system

- Clearance beneath frame for mobile hoist
- Hand switch fitted as standard on electric models
- Choice of Standard and Standard+ upholstery colours



5007 Orthopaedic Full Leg Seat Module




5006 Orthopaedic Arm Support





Available Upholstery **Colours:**



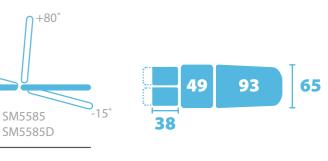

Optional	Accessories
optional	

)	5006	Orthopaedic arm support
80°	5007	Orthopaedic full leg seat module
	5008	Knee support roll set
)	5009	Orthopaedic leg support attachment
	6004	Breathing hole and plug
	6005	Paper roll holder
	6007	Removable head/knee cushion
	6052	Protective base frame cover
	6057	Battery back-up facility
	-1WF	Single foot switch for SM5565
	-2WF	Double foot switch for SM5575
	-3WF	Triple foot switch for SM5585/SM5585D




5009 Orthopaedic Leg Support Attachment



5008 Knee Support Roll Set



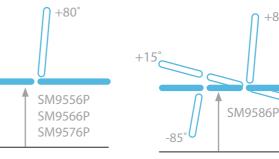

Height range 52cm to 106cm

+80

Medicare Phlebotomy Couch

Anti Bacterial 240Kg Low Transfer Height

A specialist couch designed to provide a comfortable and easy to mount seat position for blood sampling procedures. All models are fitted with a pair of multi-adjustable arm rests and side support loops for the stability and safety of the patient when mounting and dismounting the couch. The seat section is inclined for additional patient security.

Models

SM9556P SM9566P SM9576P SM9586P

Electric height, tilt and backrest



Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)
- Electric Trendelenburg tilt to -15° on SM9586P
- Patient hand-grip supports to aid mount/dismount from the couch as well, as providing patient security
- Fully adjustable arm supports for use on either side of couch
- Paper roll holder fitted as standard
- Choice of hydraulic or electric height elevation
- Excellent height range from 52cm to 106cm









Optional Accessories

6057 6007

-1WF

-2WF

-3WF

6052 Protective base frame cover Battery back-up facility Head cushion (Included on SM9586P) Single foot switch for SM9566P Double foot switch for SM9576P Triple foot switch for SM9586P

Low Patient Transfer Height





Medicare Echocardiography Couch

Anti Bacterial 240Kg Low Transfer H

Designed specifically for Echocardiography, this couch has electric operation of height, tilt and backrest angle. The backrest has a cut-away design for improved application of the transducer during scanning procedures and a specially designed stool for the clinician is also provided. The couch is also supplied with an integral footrest for the clinician, as well as a patient hand-grip and head cushion for the patients comfort and security.

Models

S	Ν	3	5	9(0	
S	М	2	6	0	h	

Electric height, tilt and backrest. Electric height, tilt and backrest. Central locking wheel system

Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)
- Electric height, tilt and backrest adjustment

- Electric trendelenburg tilt to -15°
- incorporating anti-trap design
- Ergonomic backrest designed for improved transducer positioning
- Integral footrest for clinician when traditional right-hand scanning, together with specially designed stool
- Hand-grip and side support with cushion for patient positioning
- Purpose made precision bearings on all pivot points
- Anti-bacterial powder coated frame
- Paper Roll Holder and Head Support Cushion fitted as standard
- Retractable wheel design with adjustable foot for stability on
- Choice of Standard and Standard+ upholstery colours



Available Upholstery **Colours:**






Model SM3590 with Tilt applied and backrest adjusted to ~30°



Model 6123 optional adjustable arm support for clinician



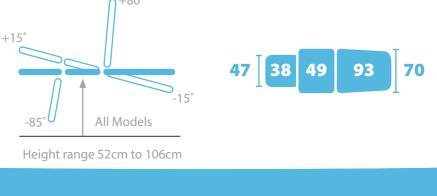

Medicare Ultrasound Couch

Anti Bacterial 240Kg Low Transfer Height

Providing smooth electric positioning of height, tilt and backrest angle, these couches feature a compact base frame and are ideal for both supine and seated procedures. The deluxe model includes the base frame cover and matching operators chair. The narrower foot section enables the fitting of optional leg supports or stirrups for gynaecological procedures.

Models

SM9686	Electric height, tilt and backrest, gas assisted
	foot section
SM9686D	Deluxe Model with Base Cover & Operators
	Chair
SM9696	Electric height, tilt, backrest and foot section
SM9696D	Deluxe Model with Base Cover & Operators
	Chair



Model 6044 optional lithotomy stirrups



Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)
- Electric adjustment of height, tilt, backrest and foot section
- Electric backrest adjustable from 0° to +80° Adjustable angle footrest from 0° to -85°
- Electric trendelenburg angle to -15°
- Choice of gas assisted or electric operated foot section
- Paper Roll Holder and Head Cushion fitted as standard
- Compact frame design with retractable wheel system
- Adjustable foot for stability on uneven surfaces
- Purpose made precision bearings on all pivot points Anti-bacterial powder coated frame
- Clearance beneath frame for mobile hoist
- 70cm standard width upholstery



MC6158 Surgeons & Sonographers Chair

Available Upholstery **Colours:**




Optional Accessories

6004	Breathing hole and plug
6041	Side support loops
6044	Lithotomy Stirrups
6045	Gynaecology Leg Supports
6052	Protective base frame cover
6054	Telescopic IV pole
6057	Battery back-up facility
-3WF	Triple foot switch for SM9686/SM9686D
-4WF	Quad foot switch for SM9696/SM9696D



Model SM9686 in trendelenburg position



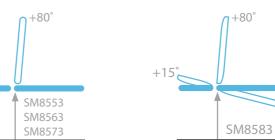

Height range 52cm to 106cm

Medicare Gynaecology Couch



Designed for general gynaecology procedures, these models are available with either hydraulic or electric height operation and feature a wide choice of optional accessories to meet both financial and clinical requirements. All models feature a crescent shaped seat design and are supplied with fully adjustable leg supports, paper roll holder, side support loops and debris tray for gynaecology procedures.

6051

Models

SM8553 SM8563 SM8573 SM8583

6043

Removable Leg Extension

Hydraulic with gas assisted backrest Electric with gas assisted backrest Electric height and backrest Electric height, tilt and backrest



Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)
- Crescent shaped seat design
- Fully adjustable patient leg supports
- Fitted with debris tray, Side Loops & Paper Roll Holder
- Choice of hydraulic or electric height elevation
- Excellent height range from 52cm to 106cm
- Adjustable foot for stability on uneven floors
- Purpose made precision bearings on all pivot points
- Anti-bacterial powder coated frame
- Electric trendelenburg tilt to -15° on model SM8583
- 65cm standard width upholstery
- Choice of Standard and Standard+ upholstery colours




ST6100 **Operators Chair**

STANDARD STANDARD+ **Available** Upholstery **Colours:** - White —— Light Grey —— Cream —— Sky Blue - Canard Black Olive ----- Lotus Green ----- Pepper Dark Blue Dark Grev






Model SM8563 Colour light grey

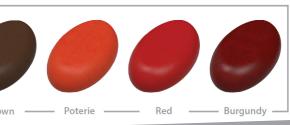



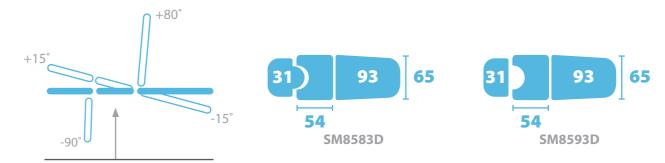

Electric control of height, backrest and tilt angle Excellent height range from 52cm to 106cm

Standard Features

Crescent shaped seat design

Fully adjustable patient leg supports

Medicare Deluxe Gynaecology Couch

Designed for colposcopy and general gynaecology procedures, these Deluxe models feature electric positioning for height, backrest angle and tilt position. They are supplied with a special set of accessories which include a matching operator's chair, base frame cover, a choice of leg extension, leg supports, matching head support cushion, side support loops, debris tray and paper roll holder.




Anti Bacterial

Transfer

240Kg

Available Upholstery **Colours:**

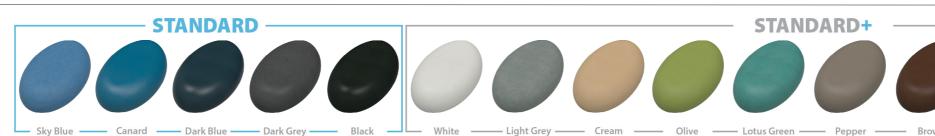









Optional Accessories





Height range 45cm to 98cm

Height range 45cm to 98cm

Medicare Bariatric 2 Section Couch

Anti Bacterial 325Kg

The Bariatric 2 Section Couches feature electric height and backrest operation for effortless patient positioning for various medical examinations and treatment procedures. There is a choice of widths available - Models SM2573 and SM2673 are available with a width of 80cm and Models SM4542 and SM4642 which have a wide frame, are available at 105cm as standard. All couches provide a safe working load and lifting capacity of 325Kg (715lbs).

Standard Features

- Safe working load and lifting capacity of 325Kg (715lb
- Choice of retractable wheels or central locking castors
- Electric backrest adjustable from -15° to +89°


- Anti-bacterial powder coated frame
- Clearance beneath frame for mobile hoist
- Hand switch fitted as standard
- Choice of Standard and Standard+ upholstery colours





Available Upholstery Colours:




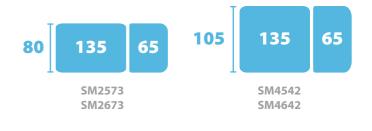

os) s	Opt 6004 6005 6074 9002P 6076 -2WF	ional Accessories Breathing hole and plug Paper roll holder Side support rail (Pair) Side support rail cushions (Pair) Memory foam upholstery (60mm) Double foot switch
5		

Model SM4642 Colour Sky Blue




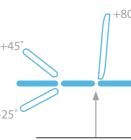

Height range 45cm to 98cm






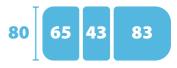







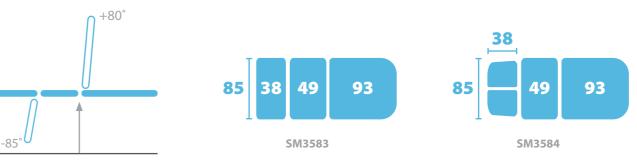

Anti Bacterial (325Kg)

The Bariatric Drop End Couches feature electric independently operated foot sections for effortless positioning for various lower limb procedures. The electric backrest supports the patient in the sitting position. There is a choice of models with either single or dual foot sections, all with a safe working load and lifting capacity of 325Kg (715lbs).

Models

SM3583Electric height, backrest and single foot
sectionSM3584Electric height, backrest and dual foot
sections

Model SM3584 Colour Dark Grey

Medicare Bariatric Drop End Couch

Standard Features

- Safe working load and lifting capacity of 325Kg (715lbs)
- Memory foam upholstery (60mm) for additional patient
- Retractable wheel design with adjustable foot for stabil
- Electric backrest adjustable from 0° to +80
- Electric foot section(s), adjustable from 0° to -85
- Electric height elevation from 50cm to 105cm
- Purpose made precision bearings on all pivot points
- Anti-bacterial powder coated frame
- Clearance beneath frame for mobile hoist
- Hand switch fitted as standard
- 85cm standard width upholstery
- Choice of Standard and Standard+ upholstery colours

Model SM3583 Colour Dark Blue.



Dual Adjustable Electric Footrests on SM3584



Available Upholstery Colours:



) t comfort lity

Optional Accessories

6005	Paper roll holder
6007	Head support cushion
6041	Side support loops (Pair)
6073	Side support rails (Pair)
9002P	Side support rail cushions (Pair)
-3WF	Triple foot switch for SM3583
-4WF	Quad foot switch for SM3584



74	160	• 74	195
_	ST7500 ST7600		ST770 ST780

Shower Trolley

Designed to provide transport to and from the showering area, these mobile Shower Trolley models feature a foam padded PVC liner for added patient comfort, a linked wheel system with central locking and steering facility, as well as fold down side support rails for patient transfer. The Junior model has an overall length of 160cm and both the Standard and Junior models are available with a rechargeable battery option as an alternative to the hydraulic height operation.

eable battery

Models

ST7700	Standard hydraulic
ST7800	Electric with Rechargeable
ST7600	Junior hydraulic
ST7500	Junior electric with Rechar

Model ST7700 Patient Transfer Height with Side & Foot Rails folded down

Standard Features

- Safe working load and lifting capacity of 240Kg (528lbs)
- 125mm central locking wheel design, with and steering facility
- Hydraulic height adjustment 41cm to 94cm via foot pedal
- Foam padded PVC liner with side supports and end panels for comfort and security
- Side support rails fold down to assist patient transfer
- Low gradient tilt towards foot end for water drainage via integral waste outlet
- Plug and drain hose supplied complete with storage bracket
- Bumpers fitted to all corners for protection during transport procedures
- Bumpers can be removed to reduce the width to 65cm to access tight doorways
- Lowers to wheelchair height for ease of patient transfer
- Clearance beneath frame for mobile hoist
- Choice of junior or standard models





Charger unit



00' 800

Height range 41cm to 94cm

- Rechargeable battery models 41cm to 94cm via hand switch
- Battery operated models include a spare battery pack and charging unit
- Purpose made precision bearings on all pivot points





Optional Accessories



-1WF Single foot switch option

-2WF

Double foot switch option (alternative to standard



-EFPS

The perimeter foot switch is a "hands free" foot operated Medicare 2 & 3 sections only.



6057

-3WF

Battery back-up available for most electric operated

6002 & 6003

Fold-down side support rails

for 2 & 3 section couches,

providing patient support


(alternative to standard triple



6052 & 6053

An easy to clean protective base frame cover that conceals the moving parts



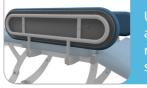

added to the side support support and protection.

If you have any queries on any of the optional accessories listed above please contact us for further details. Please note some of the accessories are supplied as standard with specific couches so please ensure you check before purchasing your chosen product as all other options come at an additional cost.



6041

Side support loops provide patient support and security and are easily removable.



The paper roll holder provides a protective barrier between the patient and the



The head support cushion is easily attached to the couch and can also be used as a support for lower limbs.



5008

The knee roll support set is an ideal accessory for supporting limbs during medical procedures.



6044

Traditional lithotomy stirrups for gynaecology procedures.





6054

Telescopic IV pole has fixing points on both sides of the couch and is easily adjusted



6004

Breathing hole & plug for improved patient comfort for prone treatments.



6006

The upper limb support has multi-angular adjustment and can be fitted to either side of the couch.



6020

A light weight free-standing adjustable leg support for lower limb procedures. Supplied with carrying case.



6045

Leg supports for gynaecology procedures. Fully adjustable for patient





Chairs & Stools



MC6158

Surgeons / Sonographers base and castors.



ST6100

Standard operators chair with 2 lever operation for height and backrest angle. Height range of 48-60cm. Black five star base and castors.



ST6103

Standard operators stool with all round release lever range of 48-60cm. Black five star base and castors.



ST6144

Economy saddle stool adjustment with a range of 48-60cm. Black five star base and castors.







ST6104

star base and glides.

MC6178

base and castors.

ST6101

Round general medical chair, standard model. Height

lever operation for height

and backrest angle. Height

foot support ring. Black five

High operators stool with all round release lever for height 53-73cm. Adjustable foot support ring. Black five star base and glides.

MC6107



and seat tilt facility. Height range of 54-71cm. Polished aluminium five star base and

Note¹: See Medical Seating brochure for more options.

Mundial Standard Range

A range of globally popular, hard-wearing vinyl colours that will retain performance for many years. Abrasion Resistance to 150,000 cycles. Available as standard on all upholstered couches and seating products.



Mundial Standard+ Range*

A selection of unique coloured, hard-wearing vinyl colours that will retain performance for many years. Abrasion Resistance to 150,000 cycles.

* Available at an additional cost for all Clinical Range Couches, excluding Shower Trolleys.



Other colours are available on request

Please Note: The illustrated vinyl colours may vary slightly from actual vinyl supplied. This occurs due to the four colour printing process, however every attempt has been made to make them as accurate as possible. 49